

Qualitative and quantitative analysis of slowly migrating chemotactic cancer cells in 3D

Nina Baumann, Verena Kuttnerberger, Elias Horn, Tina Freisinger, Roman Zantl
 ibidi GmbH, Munich, Germany, Am Klopferspitz 19, 82152 Martinsried, Germany


e-mail: tfreisinger@ibidi.de

Abstract

Chemotaxis of slowly migrating cells plays an important role in many physiological and pathological processes such as wound healing and cancer metastasis. While *in vitro* chemotaxis/migration assays of fast migrating cells such as dendritic cells are widely established, (live-cell) analysis of slowly migrating cells require a stable gradient over a relatively long period of time that is difficult to maintain.


We developed specialized chambers with confined geometries and materials that allow the investigation of slowly migrating cells along a chemotactic gradient that is stable over a period of 48h. Time-lapse phase-contrast microscopy, manual tracking and analysis of cell trajectories resulted in characteristic parameters to differentiate directed (chemotaxis) from undirected migration behavior. Our experiments revealed directed migration of HT1080 (fibrosarcoma) cancer cells in response to 10% fetal calf serum (FCS) and directed migration of MDA-MB-231 (human breast cancer) cells towards EGF in serum-free medium in a 3D Collagen gel.

A tool to study chemotaxis

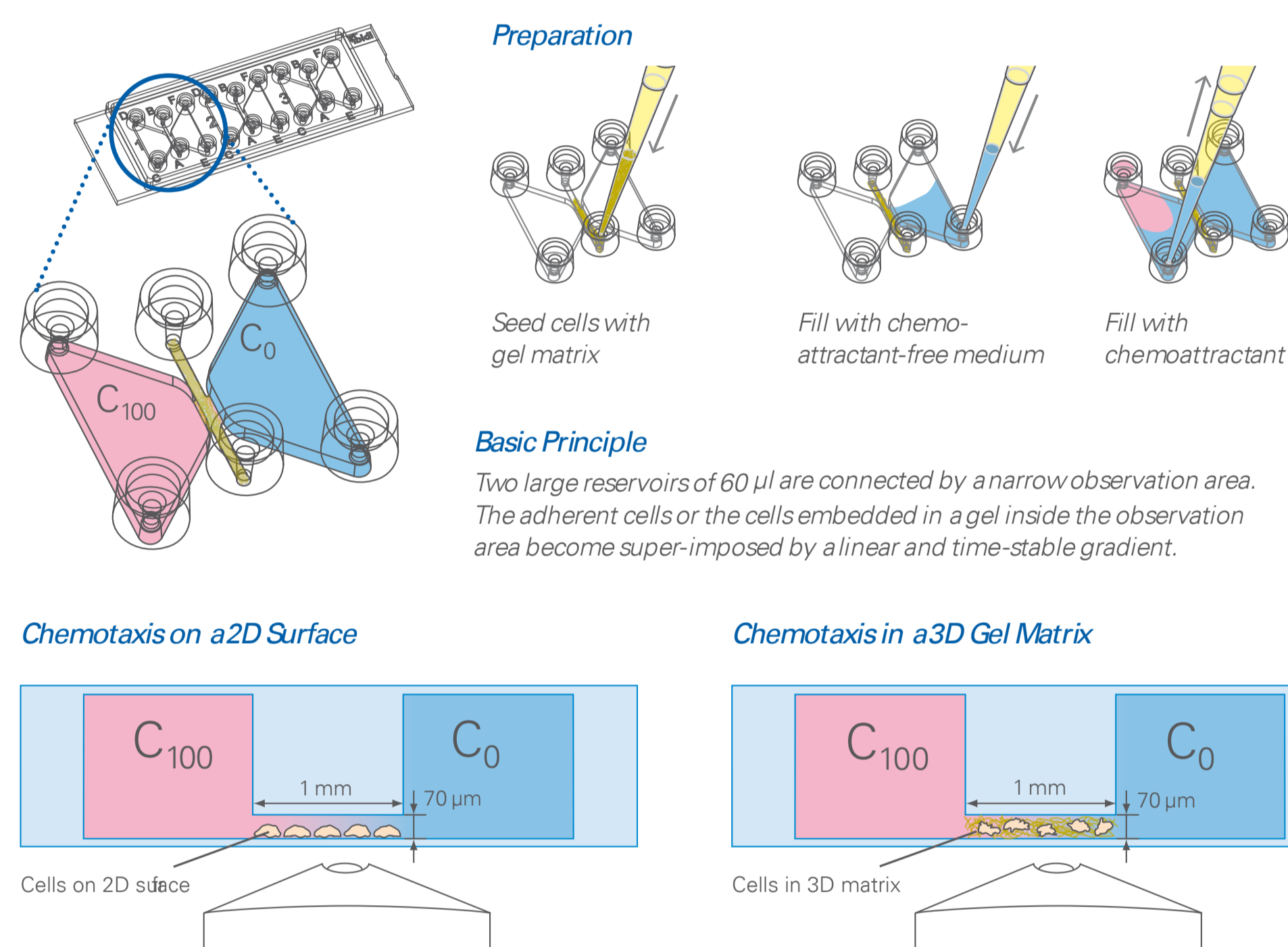


Microscopy slide to combine cell culturing and imaging
 The μ -Slide Chemotaxis is designed to analyze chemotaxis of slowly migrating cells in 2D or 3D with real-time imaging techniques. The special geometry allows formation of a chemical gradient that is stable for up to 48h.

3D chemotactic migration of HT1080 cells

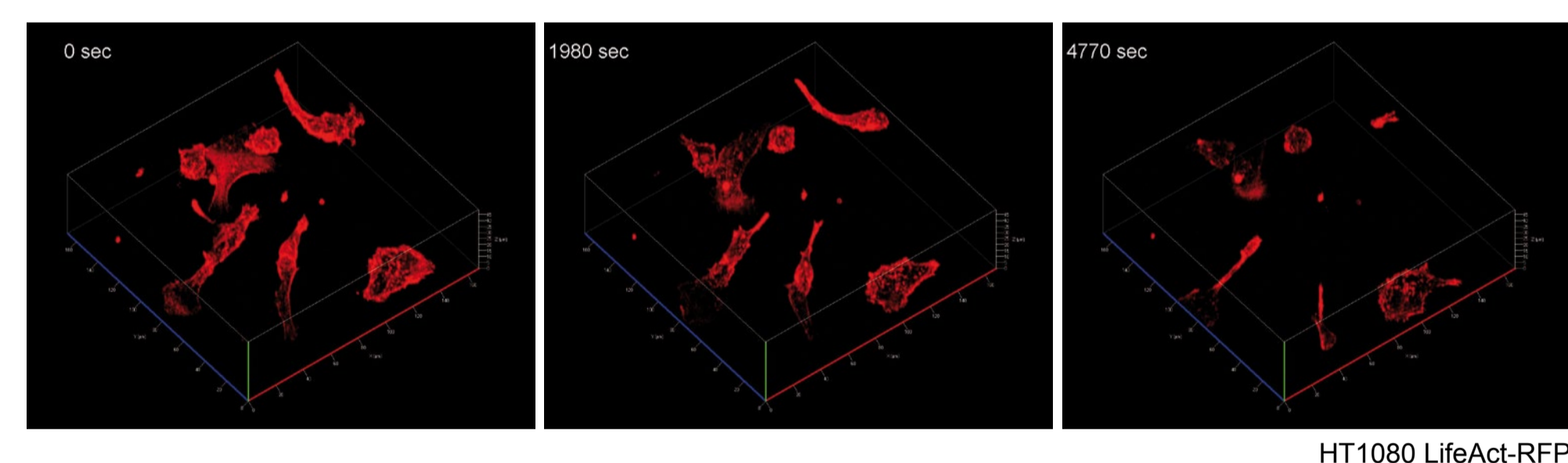

3D migration of HT1080 cells in μ -Slide Chemotaxis
 Confocal time-lapse microscopy of HT1080 cells expressing LifeAct-RFP in Collagen I gel reveals different migration behaviors dependent on substrate. Cells display ameboid-like movement on 2D surface, whereas cell shape is elongated when moving through Collagen I gel.

MDA-MB-231 cells migrate towards EGF

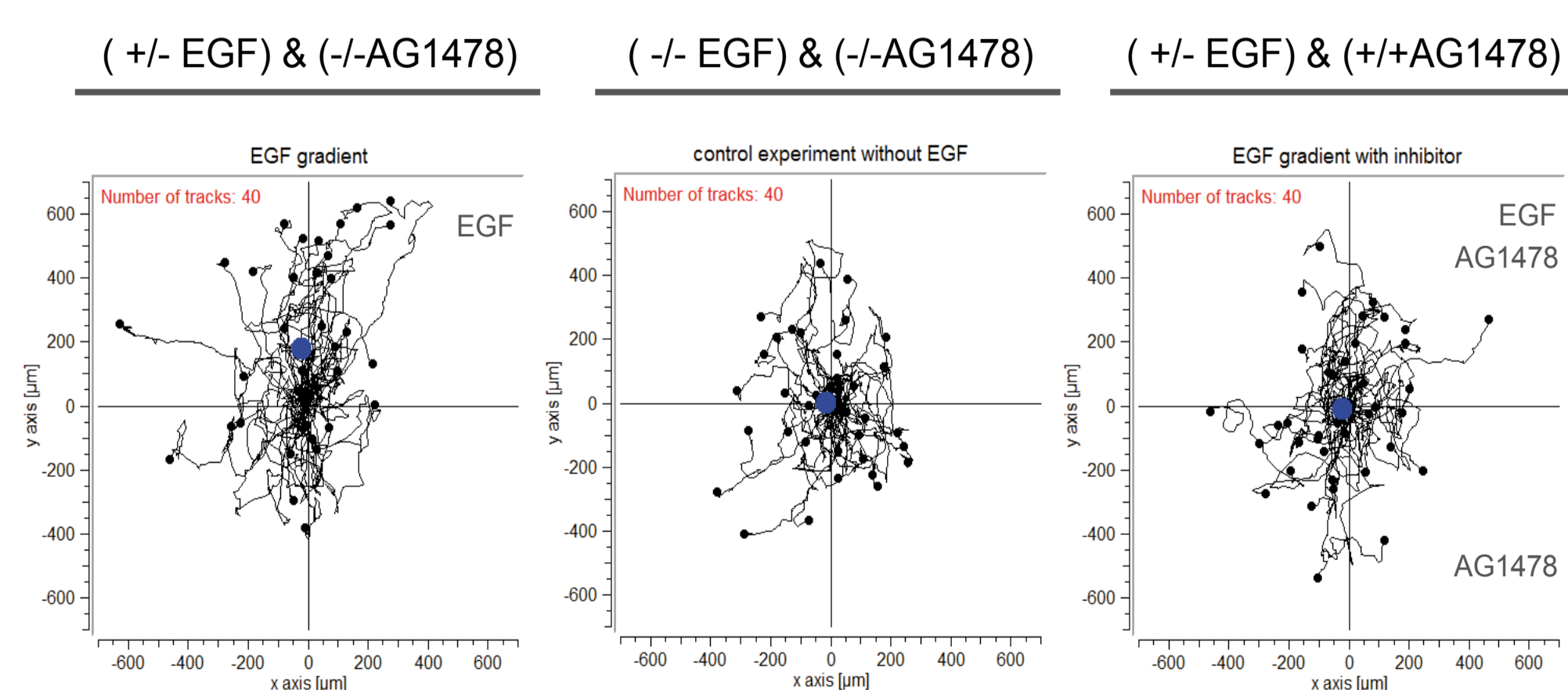

	(+/- EGF) & (-/-AG1478)	(-/- EGF) & (-/-AG1478)	(+/- EGF) & (+/+AG1478)
FMI II	0.19	0.00	-0.03
Vel	0.69 μ m/min	0.46 μ m/min	0.55 μ m/min
p	0.007	0.79	0.74

FMI ratio total distance: net distance;
 p<0.05 Distribution is inhomogeneous, hence directed

Selective directed migration towards EGF
 MDA-MB-231 breast cancer cells show directed migration behavior towards EGF in serum-free medium.

Chemotaxis of HT1080 cells towards 10% FCS

Workflow

Sample Preparation

Video Microscopy


Cell Tracking: ImageJ Manual Tracking Plugin

Visualization: ibidi Chemotaxis and Migration Tool


Quantitative and Statistical Analysis

Workflow and chemotaxis experiment of HT1080 cancer cells towards FCS
 Imaging, tracking and automated quantitative analysis of HT1080 cells reveals directed migration towards 10% FCS.


Experimental design

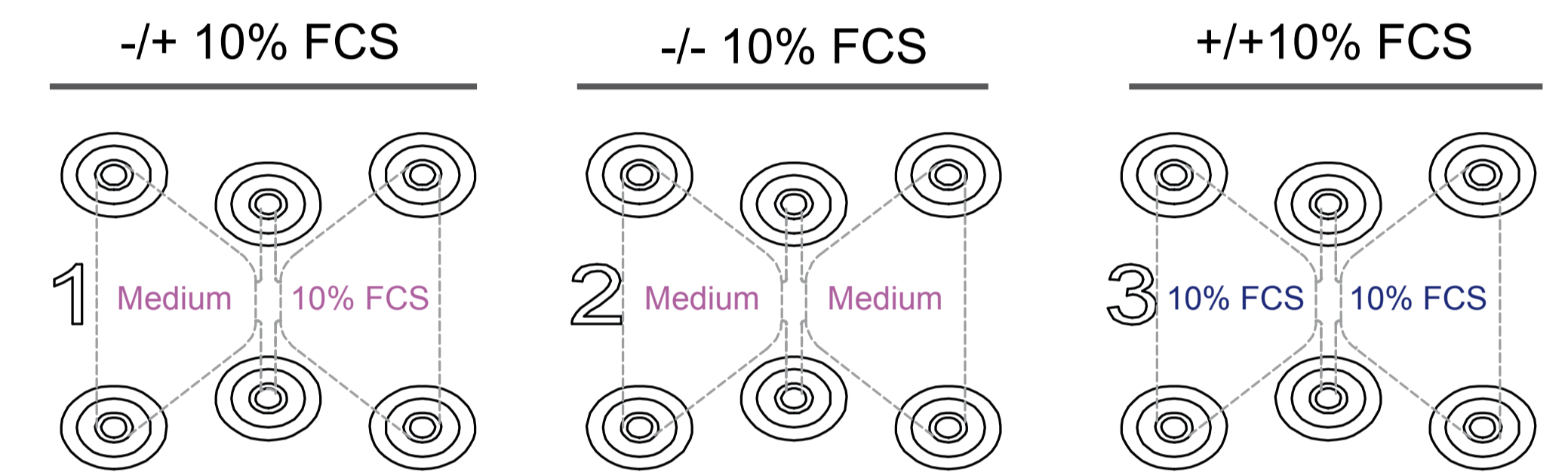


24h time-lapse, 10 min interval. 4x, phase contrast

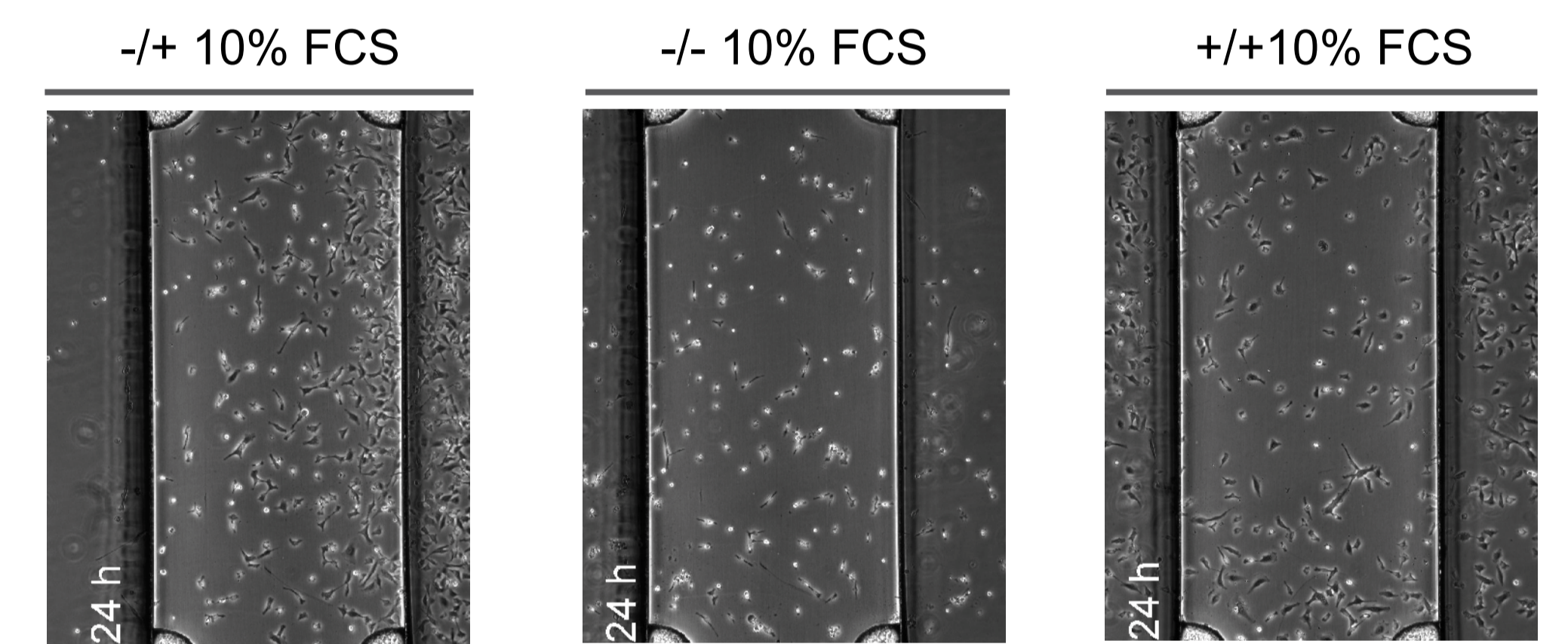

Individual frame by frame tracking of 38 cells

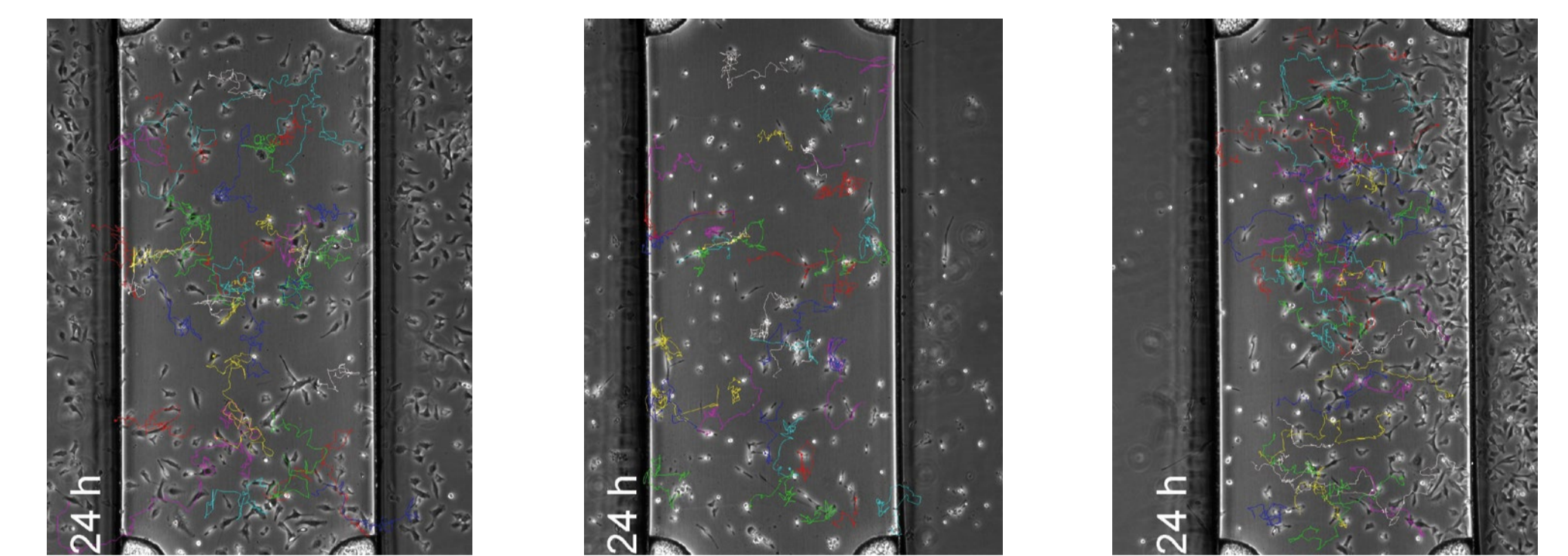

Plotting individual cell tracks and normalization to same starting point

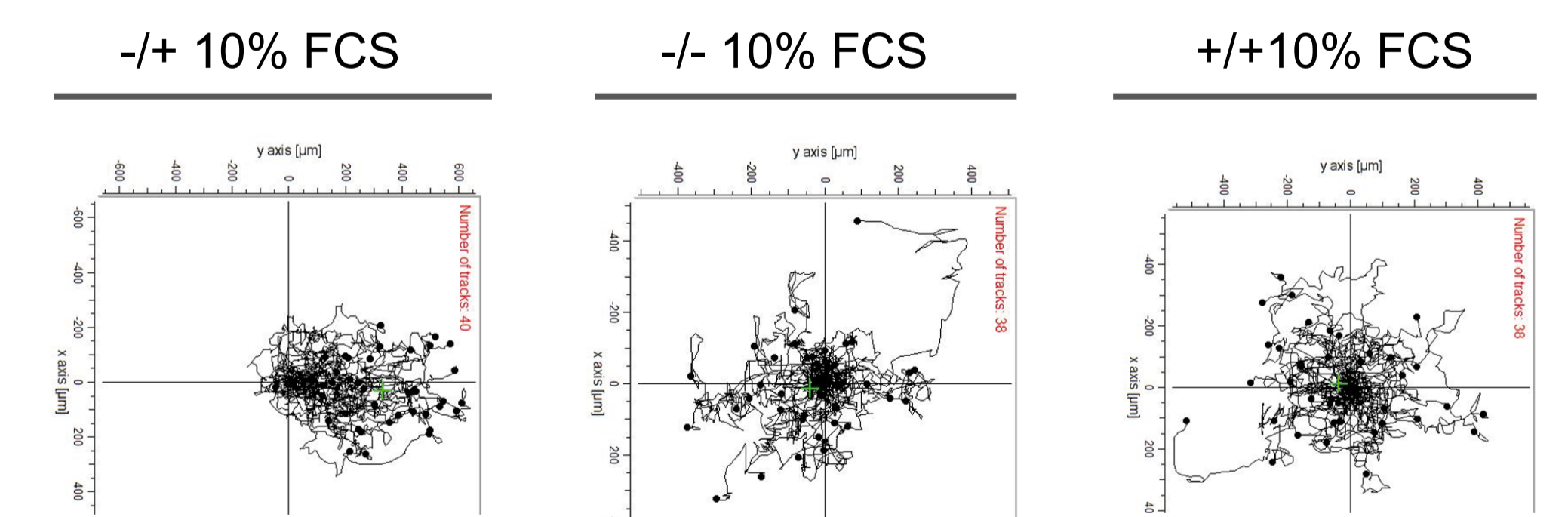

Parameter to characterize chemotactic migration

	-/+ 10% FCS	-/- 10% FCS	+/+ 10% FCS
FMI II	0.3	0.0	0.0
Vel	0.8 μ m/min	0.7 μ m/min	0.9 μ m/min
p	1.46E-13	0.351	0.302

FMI ratio total distance: net distance;
 p<0.05 Distribution is inhomogeneous, hence directed