


Part Number: **C33PU-RBD1302**

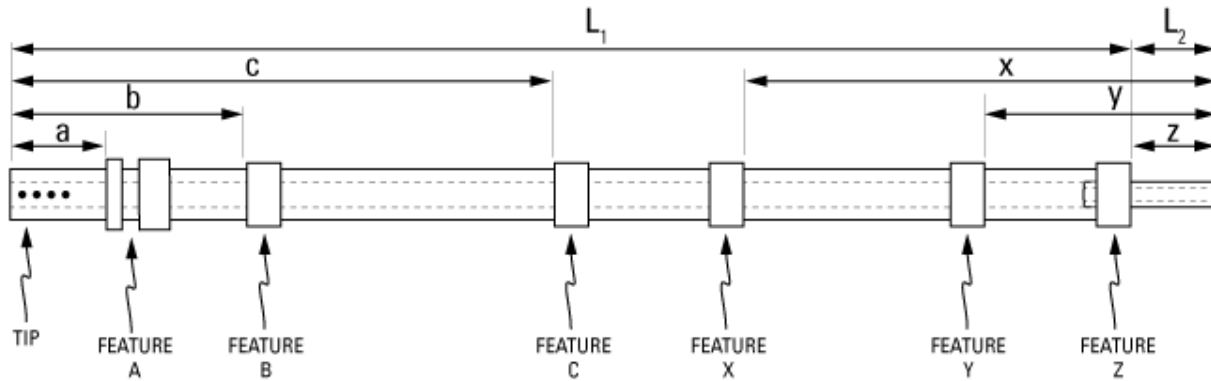

not to scale

Full Description: Catheter for rat bile duct / duodenum, 1.9-3.3Fr PU 30cm. Fits 22ga. (ReCathCo 0117EO)

Unit: pkg of 6 **Pouch Desc:** n/a **Pouch:** n/a

Approx Dead Vol: 90 μ L **Connects to:** 22ga **EO Sterilized:**

Intended Species: rat **Vessel:** bile duct

Usage Notes: catheterize duodenum with end with perfusion holes, bile duct with 1.9Fr end, cut between C and X, attach to VAHD115AB/VAHD115L bile harness

TUBING

Tip: perfusion h *Length 1* *Length 2*
 L1: 29 cm L2: 1 cm
 Material (IDxOD): PU 3.3Fr (.63x1.12 \pm .05mm) Material (IDxOD): PU 1.9Fr (.31x.64 \pm .05mm)

FEATURES

<i>Distal (tip) end</i>		<i>Location</i>	<i>Fixed</i>
Feature A:	Pair of silicone collars (2mm OD), 1.5 and 3mm wide	a: 0.8 cm	<input type="checkbox"/>
Feature B:	Large silicone collar for 3Fr (1.7mm OD), 3mm wide	b: 2.7 cm	<input type="checkbox"/>
Feature C:	Large silicone collar for 3Fr (1.7mm OD), 3mm wide	c: 12.0 cm	<input type="checkbox"/>
<i>Proximal end</i>			
Feature Z:	Large silicone collar for 3Fr (1.7mm OD), 3mm wide	z: 1.0 cm	<input type="checkbox"/>
Feature Y:	Large silicone collar for 3Fr (1.7mm OD), 3mm wide	y: 3.7 cm	<input type="checkbox"/>
Feature X:	Large silicone collar for 3Fr (1.7mm OD), 3mm wide	x: 13 cm	<input type="checkbox"/>

PARTS ADDED TO POUCH

Part 1: **Part 2:**